



Biodiesel Concentration Measurements Using Spectrum OilExpress

David Armstrong - PerkinElmer (USA), Sharon Williams - PerkinElmer (UK), Dave Wooten - Wooten Consulting
PerkinElmer 710 Bridgeport Ave. Shelton, CT 06484 Phone: 1-203-925-4602 Fax: 1-203-925-4654



Introduction

Reducing our dependence on fossil fuels and our reliance on oil and petroleum supplies are worldwide issues. Many see increasing the use of biodiesel fuel as a key initiative to meet these global needs. However, the move to include proportions of Biodiesel in everyday fuel has created a host of unresolved issues for both engine manufacturers and diesel consumers. Uppermost among these are questions concerning the concentration of the biofuel (Fatty Acid, Methyl Ester – FAME) and its quality. This application note describes how infrared transmission measurements can be used to address the concentration measurements.

Biodiesel fuels are often blended compositions of diesel fuel and esterified soy-bean oils, rapeseed oils or other potential vegetable oils, as well as fats. The physical and combustion properties of these biofuels have allowed them to achieve similar performance to diesel fuel. However, there are several characteristics (including cetane number, oxidation stability, and corrosion potential) that are of concern. These differences, especially the cetane reduction, require that adequate control of the biofuel concentration be implemented.

In addition, there are now tax incentives available in some parts of the world for the use of biodiesel. For example: in the USA this tax credit is presently in the form of a \$0.01 per FAME-% per gallon of fuel used. Therefore, the difference between 19% or 20% FAME in diesel fuel can result in a considerable tax value. A recent investigation of commercially available biofuel blends identified that 18 out of 50 splash-blended samples were not the specified 20% FAME value¹. It can be seen that there are financial justifications for an accurate biofuel concentration determination and characterization.

This work was performed using the Spectrum OilExpress system which consists of four elements:

- The PerkinElmer Spectrum™ 100 FT-IR spectrometer with high sensitivity, sampling speed and stability.
- A sealed transmission cell with zinc selenide (ZnSe) windows with a 100 μ m pathlength.
- The Molecular Spectroscopy Liquid Autosampler which provides unattended operation and rapid sample throughput of up to 50 samples per hour. The system is fitted with syringe pumps and is designed to handle samples with a wide range of viscosities, ensuring virtually no sample carryover (<0.1%).
- The PerkinElmer infrared quantitative software suite which allows analysis by various methodologies. These include Beer's law concentration calculations using Peak Height measurements and full Principal Component Regression (PCR) chemometric analysis.

The primary advantage of this system is the ability to automate the procedure from sample aspiration through report generation, including cleaning between samples. Secondly, the infrared transmission spectra carry the most information-rich data available, enabling more robust methods to be calculated.

AFNOR Method using Beer's Law

One of the few defined methods for measuring the concentration of FAME is AFNOR NF EN 14078 (July 2004) – "Liquid petroleum products - Determination of fatty acid methyl esters (FAME) in middle distillates - Infrared spectroscopy method".²

The principle of the AFNOR method is the application of a simple quantitative model of FAME content using the 1745 cm^{-1} carbonyl absorbance. When using the AFNOR methods, samples are diluted in cyclohexane to a final analysis concentration of 0 – 1.14% FAME. This produces a carbonyl peak intensity range between 0.1 – 1.1 Abs, using a 0.5 mm cell pathlength. The peak height of the carbonyl band at or near 1745 cm^{-1} is measured to a baseline drawn between 1820 and 1670 cm^{-1} . This peak height is used with a Beer's Law plot (absorbance versus concentration) to develop the calibration curve used for calculating the unknown concentrations.

While it is possible to achieve good concentration measurement, the disadvantages of this method are the need for sample dilution and the inability of the simple methodology to cope with variances in the source of the biofuel. An improved solution utilizes the more common 100 μ m flow-cell, avoiding the sample dilution errors. With the potential for increasing variance in feedstocks used to produce the FAME (namely: soybean, rapeseed or yellow-grease), peak area is proposed as a preferred calculation technique.

Peak Area Method

The modifications of this method that were employed in this study included:

- Cell pathlength – 0.1 mm
- Peak area calculation – range: 1820 – 1670 cm^{-1} with baseline set at the same range
- No Dilution – samples were not diluted to allow for the determination of the usable range

For a concentration method to be valid, the peak maximum cannot exceed the absorbance range of the spectrometer. Figures 1 and 2 demonstrate that the Beer's Law curve for this spectral region is limited to approximately 18% FAME.

In this study we took a baseline as defined in the AFNOR method at 1820 - 1670 cm^{-1} and a peak area in the same range. The sample concentration range for this method was B0 to B16 (0% to 16% FAME). The method produced a linear graph with a correlation coefficient of 0.9988. Calculating the concentration of the standards by the method yielded a Pearson's correlation of 0.9990 and a standard error of prediction (SEP) of better than 0.30%. These results indicate an acceptable method for the quantitation of FAME up to B16.

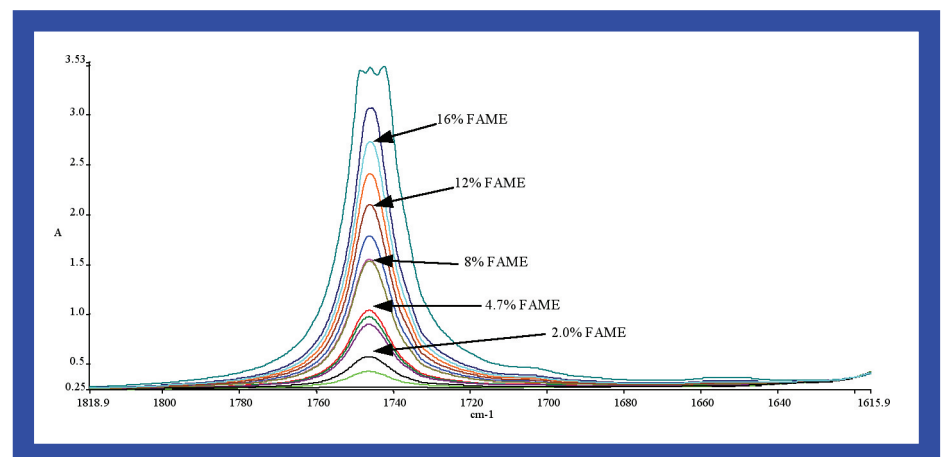


Figure 1: FTIR Spectra of varying FAME concentrations in diesel fuel.

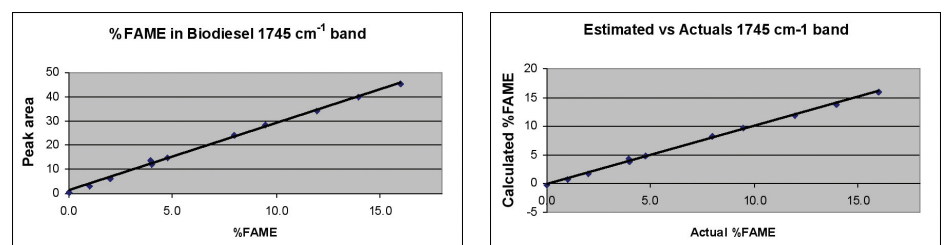


Figure 2: Beer's Law Calibration Method for 1745 cm^{-1} Peak

Further analysis of the FT-IR spectra shows additional spectral regions attributed to the FAME chemistry; for example 1300–1130 cm^{-1} (see Figure 3). The peak maximum for this spectral region does not exceed the system absorbance limit even at 49% FAME. The associated Beer's Law method uses the peak area between 1300 and 1130 cm^{-1} . Figure 4 shows the capability of this method for an extended sample concentration range from B0 to B50. The method produced a linear correlation with a correlation coefficient of 0.9997 and a standard error of prediction (SEP) of 0.38%. This is a capable method for the determination of a wider range of FAME concentrations.

Principal Component Regression Method

The peak area model is able to yield very capable calculations of the FAME concentration using short ranges of the full IR spectrum. To fully utilize all the relevant information from the whole spectrum, we moved to a chemometric analysis. In this case we used Principal Component Regression (PCR) to provide a more robust concentration assay. Samples with varying FAME

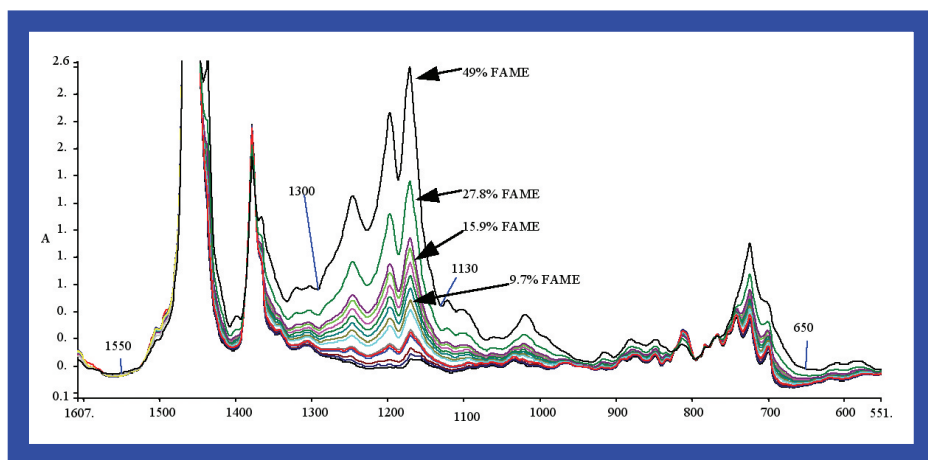


Figure 3: Fingerprint region of FAME/diesel samples.

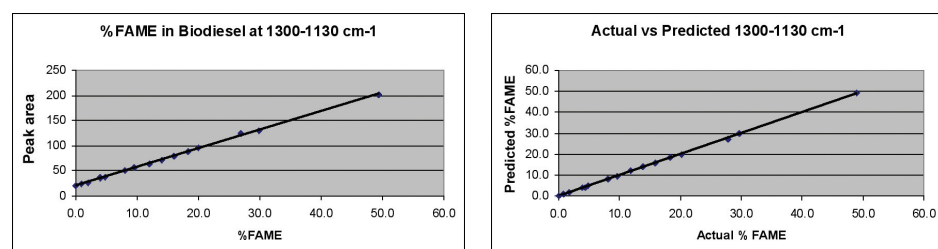


Figure 4: Beer's Law Calibration Method for 1300 – 1130 cm⁻¹.

concentrations between 0 and 20% were used in the calibration of the PCR model. The model employed as much of the entire spectrum as available. The quantitative prediction utilized only one principal component (the Regression Spectrum for the method). This spectrum (Figure 5) shows that the entire spectrum was used except the top of the 1745 cm⁻¹ FAME carbonyl peak and the C-H peaks at 2900, 1460 and 1370 cm⁻¹ region.

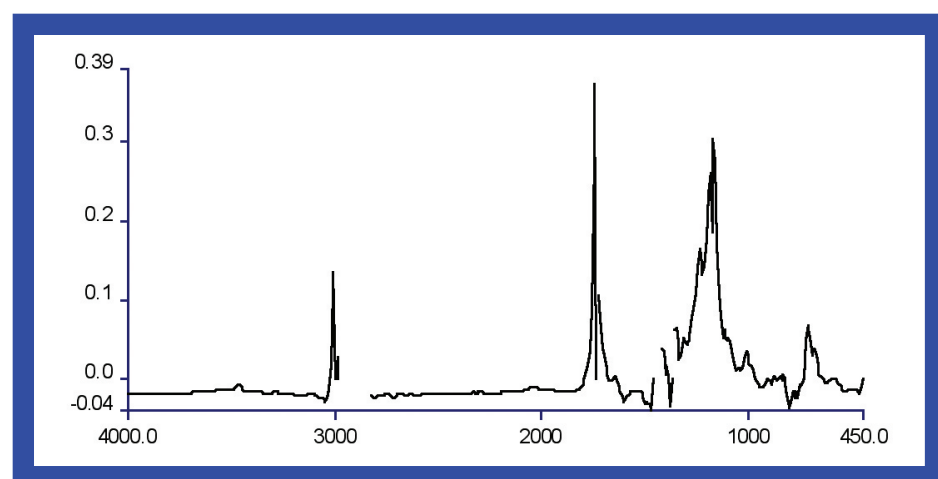


Figure 5: FAME PCR Regression Spectrum.

By using the entire spectral region, a more robust model can be generated. The statistics of this model showed a correlation coefficient of 0.9995, Pearson's correlation of 0.9997 and SEP of 0.17%. The actual against predicted results for this model is shown in Figure 6 also confirm a good prediction model.

This chemometric approach to the analyses is equal to or better than the Beer's Law methods. Although this modeling method for developing a calibration of the concentration of FAME in a biodiesel is more difficult to design, it is more robust over larger concentrations. Additionally it will allow extending the calibration range with additional samples to even higher concentrations.

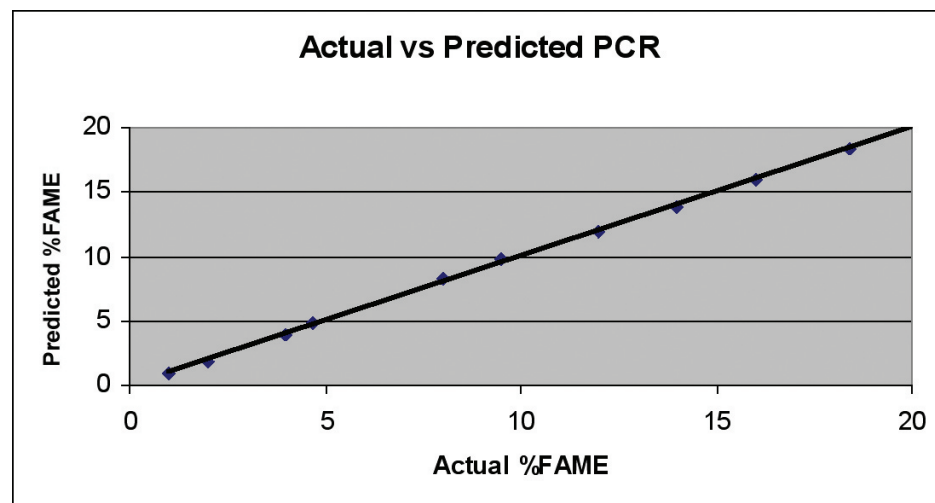


Figure 6: PCR calibration Method.

Conclusions

We have shown how infrared transmission techniques can be used to address FAME concentration measurements. All the methodologies presented achieve a standard error of prediction of less than 0.4%. This compares well with the concentration measurement of FAME in a typical "splash blend" operation, where an error of 0.5% is usually acceptable. Data analysis using either Beer's Law or Principal Component Regression (PCR) is capable of meeting this requirement.

A key advantage to using the transmission cell sampling method is that it allows auto-sampling, which can ease the routine laboratory's manpower needs. The choice of either Beer's Law or chemometrics will be determined by the particular situation. The Beer's Law approach, using peak area, benefits from being a simpler approach and is recommended for situations where there are relatively few standards and low throughput of samples. The chemometrics approach has the advantage of being more robust with respect to known constituents in the blend, better handling of interferences and reducing effect from noise contributions. Overall, PCA offers higher confidence in the quantitative prediction than is found the Beer's Law methods.

Note:

While the procedures provided in this application note may not have yet found their way into methodologies set by standard organizations or government agencies, they have been fully tested and have been demonstrated to provide quality data in numerous laboratories performing routine FAME analysis.

References:

1. NREL Technical Report TP-540-38836, Survey Of The Quality and Stability of Biodiesel and Biodiesel Blends in the United States in 2004; R.L McCormick, T.L. Alleman, M. Ratcliff, L. Moens, and R. Lawrence October, 2005
2. AFNOR NF EN 14078 Liquid Petroleum Products - Determination of Fatty Acid Methyl Esters (FAME) in Middle Distillates - Infrared Spectroscopy Method, July 2004.
3. ASTM Proposed Standard Test Method for "Determination of Biodiesel (Fatty Acid Methyl Esters) Content in Diesel Fuel Oil Using Mid Infrared Spectroscopy"

Acknowledgements:

This work was performed for PerkinElmer by David L. Wooton, PhD, of Wooton-Consulting, Beaverdam, VA, USA.

New Analyser for In-the-field Biodiesel Screening

Zeltex (USA) is proud to announce its ZX-101XL portable fuel analyser is now able to screen biodiesel in the field and at the pump. For over thirteen years, state and federal governments, as well as numerous private companies have been using the Zeltex ZX-101 line of fuel analysers to accurately test octane and cetane in the field. Now the ZX-101XL can perform the same test on biodiesel fuel. With calibrations for biodiesel percentage and ethanol percentage, the ZX-101XL will prove to be the only choice for in-the-field fuel screening. Across the United States and in forty countries worldwide, Zeltex's fuel analysers have established themselves as the analyser of choice. Their analysers will provide you with lab-accurate and dependable readings within sixty seconds. Operating on "AA" batteries, the ZX-101XL can be used to test biodiesel and ethanol-blended fuels.

New, Specifically Designed Instrument to Analyse the Oxidative Stability of Biodiesel, Diesel And Blends

For the last year **Metrohm** (Switzerland) have been developing a new instrument for analysis of diesel and biodiesel. They have improved on the existing technology of the 743 Rancimat to specifically design an instrument for the biofuels market.

During the measurement a stream of air is passed through the sample which is contained within a sealed and heated reaction vessel. These conditions accelerate the rate of oxidation of the fatty acid methyl esters in the sample, with peroxides being initially formed as the primary oxidation products. After some time the fatty acid methyl esters disintegrate completely; the second oxidation products formed include low-molecular organic acids in addition to other volatile organic compounds. These are transported in a stream of air to a second vessel containing distilled water. The conductivity in this vessel is recorded continuously. The organic acids can be detected by the increase in conductivity. The time elapsed until these secondary reaction products appear is known as the induction time or induction period.

The 873 has been modified to manage the very aggressive nature of biodiesel. Modifications include chemically resistant iso-versinic tubing and chemically resistant glass measuring vessels.

Metrohm have also proved that this principle method can be used to measure the oxidative stability of diesel and biodiesel/diesel blends. When developing the method, it was found that with volatile oils, mineral diesel and diesel blends there was a degree of evaporation on heating, resulting in wrong determination values. This was overcome by developing longer (250 mm) reaction vessels for measurements with these particular samples. These longer vessels reduce evaporation loss through the refluxer condensation principle which occurs within the longer vessels giving accurate reproducible results.

The updated software includes flexible and comfortable user administration with login functions; this can be used to define detailed access rights for groups of users and individuals. The software parameters are now also specific for biodiesel and diesel analysis.